


# A Mobile Aneroid That Measures Up To Our Standards

## The Roll-By® Mobile Aneroid from W.A. Baum

The Roll-By® Aneroid is built to agree with our master instruments, which are traceable to the United States National Institute of Standards and Technology. It complies with ANSI/AAMI SP-10 for accuracy and performance (1% full scale  $\pm$  3 mmHg)

*Heavy gauge, machined aluminum tubing*

*Designed to roll quietly, wheels are 2" (50mm) hooded casters with mar-resistant Polyurethane treads*

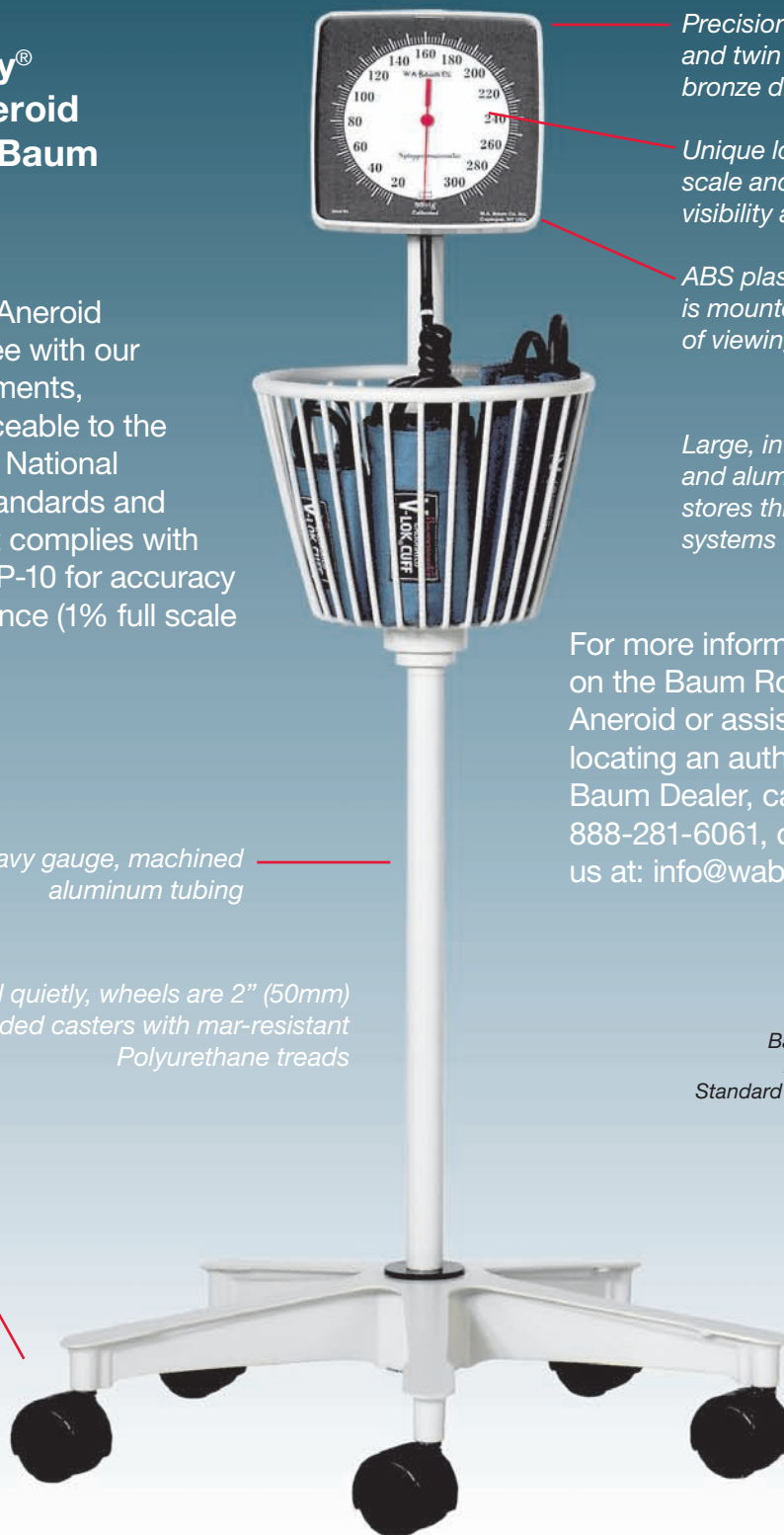

*Precision low-friction movement and twin pre-seasoned phosphor-bronze diaphragm capsules*

*Unique low reflection, high contrast scale and red pointer increase visibility and reduce parallax error*

*ABS plastic manometer housing is mounted at a 20° angle for ease of viewing.*

*Large, integral, nylon coated steel and aluminum storage basket easily stores three complete inflation systems (available latex-free)*

For more information on the Baum Roll-By® Aneroid or assistance in locating an authorized Baum Dealer, call toll-free 888-281-6061, or contact us at: [info@wabaum.com](mailto:info@wabaum.com)

*Baumanometer®.  
Blood Pressure  
Standard the World Over*



## W.A. Baum Co. Inc.

Since 1916

Originator and Maker of Accurate Blood Pressure Apparatus Exclusively

620 Oak Street, Copiague, NY 11726 U.S.A.

[www.wabaum.com](http://www.wabaum.com) • [info@wabaum.com](mailto:info@wabaum.com) • Toll-Free: 1-888-281-6061 • Tel: 631-226-3940 • Fax: 631-226-3969

# Features and Specifications



1820

## Calibrated® V-Lok® Inflation System:

**Consists of:** Bulb, Air-Flo® Control Valve, Inflation Bag and Calibrated® V-Lok® Cuff.

**Accurate:** The Calibrated® V-Lok® System immediately indicates whether the correct size cuff is being used. The center of the bladder is clearly marked on the inner and outer surface of the cuff for correct application on either the right or left limb.

**Comfortable:** Cuffs conform to the shape of the limb and have no hard, stiff edges.

**Durable:** Baum cuffs are made of tightly woven, heat set Dacron™ polyester fabric. They are double stitched at critical points for long wear and treated with an effective anti-microbial agent. Hook and loop fasteners of matched holding strength provide the optimal number of open-close cycles.

**Non-Latex:** Available for those who are hypersensitive to natural rubber latex.

**Guaranteed:** Latex and **Non-Latex** inflation bags, bulbs and tubing are guaranteed for five years.

For additional information, refer to the Calibrated® V-Lok® Cuff Data Sheet or visit [www.wabaum.com](http://www.wabaum.com).

## Roll-By® Mobile Aneroid Description:

- 300 mmHg clinical sphygmomanometer for mobile use.
- Complies with ANSI/AAMI SP-10 for accuracy and performance (1% full scale ± 3 mmHg).
- Available with Latex or **Non-Latex** components.
- Unique, low reflection, high contrast scale and red pointer to increase visibility and reduce parallax error.
- Diecast aluminum 5-wheel base (20" dia.) with 2" hooded casters and 2 lb. cast iron weight for increased stability.
- Manometer housing is mounted at a 20° angle for ease of viewing.
- Heavy gauge machined aluminum tubing finished in Neutral Gray baked enamel.
- **Calibrated® V-Lok® Inflation System** and 8 feet of Extendex® coiled tubing.

## Baum Roll-By® Model 1150 - Large Basket Specifications:

The **Baum Roll-By® Model 1150** has a large circular basket to store a full complement of inflation systems.

**Instrument Finish:** Neutral White ABS plastic

**Standard Cuff Color:** Medium Blue

**Packaged Dimensions:** 22" L x 21½" W x 15¼" H

**Packaged Weight:** 17.70 lbs.

## Roll-By® Model 1150 with 8' Extendex® tubing and large circular basket.:

Cat. No.	Non-Latex	Cuff Size(s)	Inflation System	Cuff Holder
► 1150	1150NL	Adult	Calibrated® V-Lok®	Large Circular Basket
1151	1151NL	Child/Small Adult	Calibrated® V-Lok®	Large Circular Basket
1152	1152NL	Infant	Calibrated® V-Lok®	Large Circular Basket
1153	1153NL	Newborn	Calibrated® V-Lok®	Large Circular Basket
1154	1154NL	Thigh	Calibrated® V-Lok®	Large Circular Basket
1155	1155NL	Large Arm	Calibrated® V-Lok®	Large Circular Basket
► 1170	1170NL	Adult, Child, and Large Arm	Calibrated® V-Lok®	Large Circular Basket
1140	1140NL	Roll-By® Manometer Only		

## Guarantee:

Every Baum Aneroid instrument is guaranteed against defects in workmanship and material for a period of one year from date of purchase. This guarantee is invalid when other than genuine W.A. Baum parts are used. Latex and latex-free components are guaranteed for 5 years.

We recommend that every aneroid instrument be checked for accuracy periodically by following our "Testing Instrument Accuracy" procedure\* to determine if it needs adjustment. Adjustments will be made to Baum aneroid instruments without charge when returned, postage prepaid, unless the instrument has been damaged in shipment, tampered with or abused. Our pledge is that our customers will be completely satisfied with our products.

\*Printed copies of "Testing Instrument Accuracy" are available on request or by visiting the FAQ section of our website at: [www.wabaum.com](http://www.wabaum.com).

Baumanometer, Standby, V-Lok, Extendex, Calibrated, Hi/Lo, Air-Flo and Roll-By are registered trademarks of W.A. Baum Co., Inc.

**W.A. Baum Co., Inc.**

620 Oak Street  
Copliague, NY 11726 USA

Tel: 631-226-3940

Fax: 631-226-3969

Toll Free: 1-888-281-6061

Web: [www.wabaum.com](http://www.wabaum.com)

Email: [info@wabaum.com](mailto:info@wabaum.com)

**Since 1916 Originator and Maker of Accurate Blood Pressure Apparatus Exclusively**

**W.A. Baum Co. Inc.**



1150